

# Percuplasty™ SuperBite™

Fully-Threaded Compression Screws

## Enhanced compression and stability with every turn.



✓ **Advanced compression and stability**

Powered by a continuously variable thread pitch and optimized head design

✓ **Rapid screw insertion**

Self-drilling and self-countersinking design maximizes OR efficiency

✓ **Engineered for minimal soft tissue irritation**


Beveled head\* and rounded edges provide reduced prominence

# Percuplasty™ SuperBite™

Fully-Threaded Compression Screws



Self-drilling tip

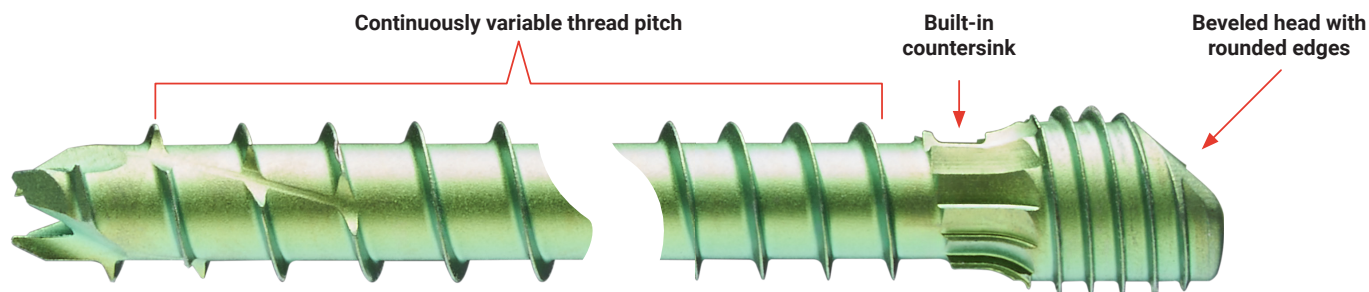

Sterile-packed implants for **convenient delivery**

	Part# (XX=Length)	Screw Diameter	Head Diameter	Lengths Offered	Wire Size
	SC25XX	2.5mm	3.0mm	10-40mm <i>in 2mm increments</i>	1.1mm
	SC35XX	3.5mm	4.2mm	12-50mm <i>in 2mm increments</i>	1.4mm
	SC40XX	4.0mm	4.8mm	24-60mm <i>in 2mm increments</i>	1.6mm
	SC50XX	5.0mm	6.0mm	30-90mm <i>in 2mm increments up to 50mm, 5mm beyond 50mm</i>	2.0mm
	SC70XX(X)	7.0mm	8.0mm	40-130mm <i>in 5mm increments</i>	3.0mm

\*Beveled head available in 3.5, 4.0, and 5.0mm screw diameters only.

Before use of the system, the surgeon should refer to the appropriate instructions for use and surgical technique for complete warnings, precautions, indications, contraindications, and adverse events. Risks include, but are not limited to: infection, pain, discomfort from the presence of the implant, loosening of the implant, and loss of correction with nonunion or malunion. If any of these occur, additional treatments may be needed. Additional information about risks, warnings, and instructions is available at [Lapiplasty.com/surgeons/labeling](http://Lapiplasty.com/surgeons/labeling). Refer to TMC Screw Fixation IFU (LBL 1405-9056) for additional information.

For more information, go to  
**Treace.com**



**TREACE**  
A Step Ahead.™